

## ***Why is PAP Clean better than cheaper models on Amazon?***

### **The Ozone Generator**

Most cheaper models only produce ozone

- Ours uses Ozone and a Room-Air Flush

Most cheaper models have only one mode

- Ours has a 30 or 20 minute cleaning cycle

Most cheaper ones don't have a visual screen

- Ours has an advanced color screen
- Ours displays the battery level
- Our displays cleaning cycle images
- Ours displays the remaining cycle time

### **The Sanitizing Bag**

Most cheaper models don't include a bag

- Ours includes the strongest bag available

Most with bags are too small for Full Masks

- Ours can hold 2 masks and humidifier tubs

Most with bags need them replaced often

- Ours does not require new bags or filters

### **Power Options**

Most cheaper models don't have a battery

- Ours has a 2,600mAh Lithium-ion battery

Most with batteries need daily recharging

- Ours runs 6+ cleanings before a recharge

Most with batteries don't include a charger

- Our includes a complete 110V power supply

## ***In life, you get what you pay for...***

Ask your CPAP professional for some feedback, patient success stories and why they recommend **PAP Clean**.

**PAP Clean**

[www.PAPclean.com](http://www.PAPclean.com)

**PAP Clean**

**Portable CPAP Sanitizer System**  
Eliminates 99% of Germs & Bacteria



### **Most Features of any CPAP Cleaner**

- Trusted Ozone Technology
- No need to replace Filters or Bags
- Advanced color screen displays each step of the cleaning process





Easily Cleans Your Mask,  
Humidifier Tank & Tubing System.  
**Works with All CPAP Brands!**



## COMPLETE SYSTEM INCLUDES:

### High Efficiency Ozone Generator

- 30 minute Ozone Cleaning cycle, then a 5 min Room-Air Flush reverts ozone, and a 60 minute resting countdown
- Screen displays each step in the process
  - Icon for Cleaning, Flush and Resting
  - Timer counts down minutes of each cycle
  -  icon when cleaning,  when complete
- Dual Sanitizing: Standard or Quick Cleaning
  - Sanitizing Mode displayed on screen
- Built-in 2,600mAh Lithium-ion Battery
  - Battery level displayed on screen

### Large Zippered Sanitizing Bag

- Doubles as a travel bag for your PAP Clean and CPAP accessories
- Can simultaneously clean multiple masks & humidifier tanks
- Tube Guard™ strap prevents ozone leaks

### No Additional Accessories Needed


- Complete system cleans ALL Heated and Travel CPAP Tubes

### Type C USB Cord and Power Supply

- Runs 6+ sanitizing cycles per charge

# Sanitize in Three Easy Steps!

1

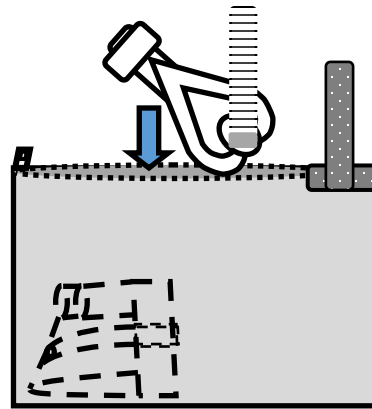


Place your CPAP Mask and empty water tub in the PAP Clean Bag.

Zip the bag closed and connect the CPAP Tube to the Ozone Generator.


Turn the PAP Clean on and walk away. The screen displays each step of the cleaning process.

2

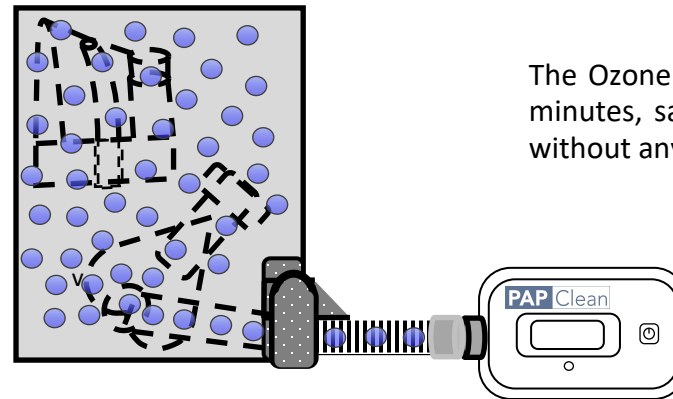

The Ozone Generator quietly runs for 30 minutes, sanitizing everything in the bag without any detergents or liquids.

3

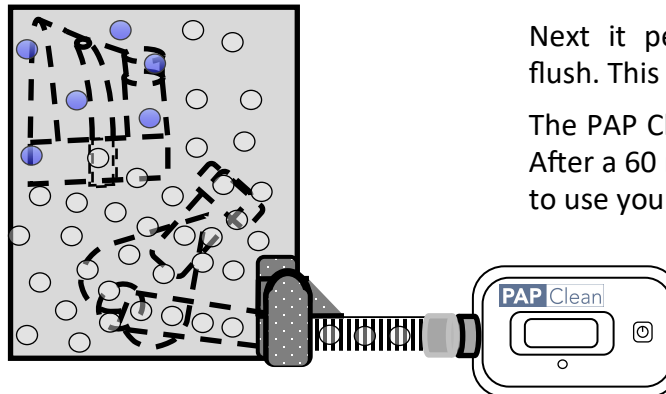

Next it performs a 5 minute room-air flush. This assists in the ozone reversion.

The PAP Clean will automatically turn off. After a 60 minute resting period, it is safe to use your CPAP.